



Yes, you can.[®]

Invacare[®] Propad[®] Mattress Overlay

The Original Mattress Overlay


Widely used by healthcare professionals in Hospitals, Nursing Homes and in the Community, the **Invacare[®] Propad[®] Mattress Overlay** has proved itself to be an essential element of a successful pressure care programme. Designed to sit on top of a conventional mattress, the **Propad[®] Mattress Overlay** provides comfort and pressure area care to those considered at 'High Risk' (Waterlow score ≥ 20)* of developing pressure ulcers.



3 year warranty

➤ Effective Pressure Ulcer Prevention

The foam core has independent, geometrically-cut surface cells which significantly reduce shear and friction, whilst increasing the overall support surface area, thereby reducing peak pressure.

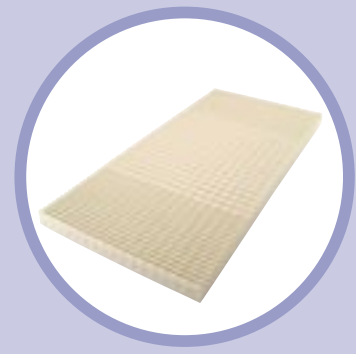

➤ Even Weight Distribution


The **Propad® Overlay** has larger, more supportive cells at the centre and smaller cells at each end. This enables the mattress overlay to achieve even weight distribution by conforming to the shape and form of the body and cradle bony prominences.



➤ Essential Air Circulation


The **Propad® Overlay** has a network of 'keyhole' channels at the base of each cell. This ensures air circulation within the overlay, minimising the build up of heat and moisture, helping to reduce shear and friction, and prevent sweating.

Features and Options

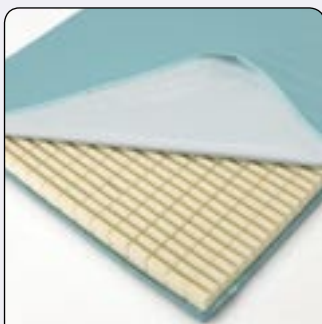


Cover helps to prevent a build up of shear and friction forces

Cover is two-way stretch, water-resistant and vapour-permeable.

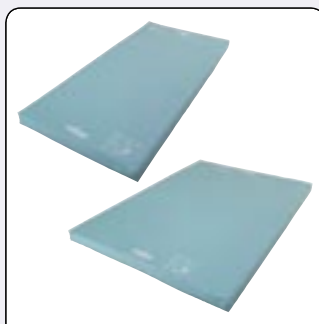

Available in a range of sizes

Single - 187.5 x 87.6 x 8.5 cm,
Double - 187.5 x 137.2 x 8.5 cm,
King Size - 198.0 x 152.4 x 8.5 cm.

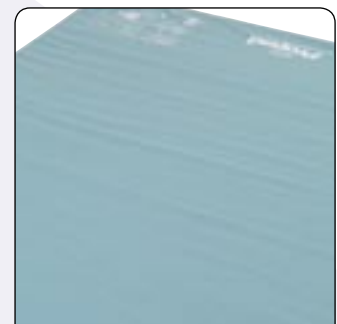


Covers have an anti-slip base


Ensures safe use of the **Propad® Overlay** and helps to prevent slippage.

Warranty
3 Years

Technical data

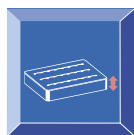

Height: 8.5 cm

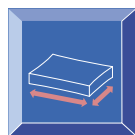


Std sizes:
187.5 x 87.6 x 8.5 cm*1
187.5 x 137.2 x 8.5 cm*2
198 x 152.4 x 8.5 cm*3

Additional sizes available on request



Max Client
Weight:
108 kg
(17 stone)

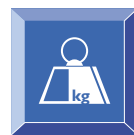

Weight of product:
5.5 kg



Cleaning: 80°C
Refer to instructions
for use and relevant
infection control
guidelines

Water Resistant
Vapour Permeable

*1 Single
*2 Double
*3 King Size

Manufactured to European Standards: BS 7177 Crib 5

Cover Colour



Green (PU)

* Essential nursing care is pivotal in pressure ulcer prevention. This Overlay will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. The Waterlow score is a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

All information quoted is believed to be correct at time of print. Invacare® reserves the right to alter product specifications without prior consultation.



INVACARE® Limited
Pencoed Technology Park · Pencoed · Bridgend CF35 5HZ
Tel: +44 1 656 776 222 · Fax: +44 1656 776 220
E-mail: uk@invacare.com · www.invacare.co.uk
Sales Order E-mail: ordersuk@invacare.com

INVACARE® Ireland Ltd
Unit 5 · Seatown Business Campus
Seatown Road · SWORDS · County Dublin
Tel. +353 1 8107084 · Fax +353 1 8107085
Email: Ireland@invacare.com · www.invacare.ie